
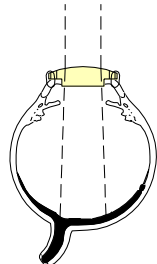


Ocular Landers Wide Field Temporary Keratoprosthesis

	Product Code	Image Mag	Contact Diameter	Static FOV	<i>Designed with: Maurice B. Landers, III, M.D. Chapel Hill, NC</i>	
	OLTK - 7.2	2.29x	7.2mm	28°		
	OLTK - 8.2	2.29x	8.2mm	30°		

Design

- The Landers Wide Field Temporary Keratoprosthesis is designed with a convex anterior surface to facilitate viewing of the peripheral retina and the posterior pole during vitrectomy and corneal transplant surgery.
- This lens has six suture holes around its peripheral edge. Lens design is suitable for phakic, aphakic and pseudophakic eye.
- Once the sutures are in place, the lens/cornea interface acts as a seal and provides a closed system for a vitrectomy procedure.
- May be used in conjunction with our handled vitrectomy lenses (OBVI, OFVI, OMVI, OPGVI, OPVI-3) or the Machemer Flat Vitrectomy Lens (OLV-5).

Technique

- A 7.0mm or 8.0mm corneal trephine is used to remove corneal button. The 7.0mm is most commonly used.
- The Keratoprosthesis (0.2mm larger than the button) is placed through the opening and rotated down into the cornea.
- The Keratoprosthesis is sutured to the cornea using 8.0 nylon suture.
- A Flieringa ring should also be sutured to eyes with flexible scleras.
- The large aperture provides a wide angle view of the posterior retina and permits clear visualization of the peripheral retina when combined with scleral depression.
- In some cases, the fundus view can be improved by placing a Machemer Flat Vitrectomy Lens on the keratoprosthesis with a drop of visco-elastic material between the lenses.
- At the conclusion of the pars plana vitrectomy, the sutures holding the keratoprosthesis are removed. The donor button is sutured in place.

Cleaning & Disinfection

See Cleaning Method 8



2255 116th Ave NE, Bellevue, Washington 98004-3039 USA
T: 425-455-5200 or 800-888-6616 F: 425-462-6669
E: contact@ocularinc.com I: www.ocularinc.com

© 2001 Ocular Instruments
Revised on 2021-06-07 5516J11435