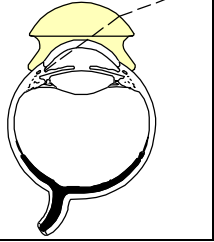



Ocular Koeppe Diagnostic Lenses

| CE | Product Code | Style | Image Mag | Contact OD | Static Gonio FOV |  |
|---|--------------|--------|-----------|------------|------------------|---|
|  | OKL | Large | 1.50x | 19mm | 160° | |
| | OKM | Medium | 1.53x | 18mm | 160° | |
| | OKS | Small | 1.57x | 17mm | 160° | |

Design

- § The Koeppe Diagnostic Lenses use direct gonioscopy to examine the angle.
- § The lens rests on the scleral band creating a corneal vault.
- § Gonioscopic solution must be employed.
- § A light source, head loupe or microscope is required to achieve desirable resolution and illumination.
- § It is necessary to recline the patient so the practitioner can view 360° of the angle by walking around the patient.

Cleaning & Disinfection

See Cleaning Method 1



2255 116th Ave NE, Bellevue, Washington 98004-3039 USA
 T: 425-455-5200 or 800-888-6616 F: 425-462-6669
 E: ocular@ocular-instruments.com I: www.ocular-instruments.com

© 2001 Ocular Instruments
 5510G3275