



OCULAR

Latina 1X Indexing PSLT Lens



Introducing the **NEW** Ocular Latina 1X Indexing Pattern Scanning Laser Trabeculoplasty (PSLT) Gonioscopy lens *featuring our revolutionary integrated indexing ring*. Provides a rotational reference for each pattern of laser spots. The lens *“clicks”* into 32 positions that coincide with the embedded reference bar within the lens. Anterior ring allows trouble-free, one handed rotation of the lens that provides a quick and precise location at each 11.25° of rotations. Can be used with all Selective Laser Trabeculoplasty procedures. The anterior ring can be easily removed for cleaning.

U.S. Patent #7,766,480

U.S. Patent #8,861,061

Product Code: **OL1PLTF**
Gonio Mag: **1.0x**
Gonio Laser Spot Mag: **1.0x**
Contact Diameter: **18mm**
Lens Height: **34mm**
FOV: **130°**


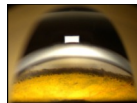


Photo represents image of a model eye

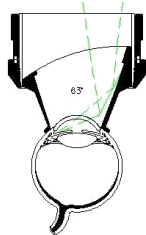

Retail: **\$975.00**

For more information on this product, including a **video demonstration on how the Indexing Ring operates**, contact our customer service department and ask for **video #11824** or visit our website at www.ocularinc.com.

Ocular Instruments, Inc.
2255 116th Ave N.E.
Bellevue, WA 98004
T: 800-888-6616 | F: 425-462-6669
contact@ocularinc.com | www.ocularinc.com



See beyond the surface