




Creating Solutions to Save Sight

Ocular Brasse Wide Field Capsulotomy Lens


For **60 years**, Ocular Instruments has been Creating Solutions to Save Sight, and is now proud to introduce the latest addition to our line of Capsulotomy YAG lenses, the **NEW Ocular Brasse Wide Field Capsulotomy Lens**. While traditional capsulotomy YAG lenses have been the standard for decades, the New Ocular Brasse Wide Field Capsulotomy lens provides increased FOV necessary to create the larger openings, required for the premium IOL and reduce dysphotopsias.

The Ocular Brasse Wide Field Capsulotomy Lens allows viewing of the anterior chamber to the anterior fundus (typical of capsulotomy lens). Larger diameter optic allows more peripheral view compared to traditional capsulotomy lens designs. Made of glass and aluminum, it is compatible with **standard steam sterilization** cycles.


Prod Code: **OBWAC**
Contact Diam: **16mm**

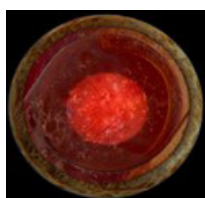

Actual view of the capsular opening using a standard capsulotomy lens

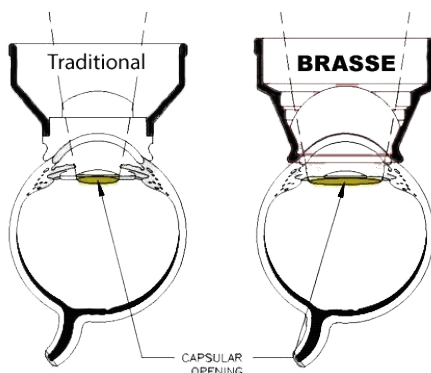

The large FOV of the Brasse Capsulotomy lens allows a wider area of the capsular membrane to be viewed and treated by the YAG laser. The larger treated area helps to reduce annoying light artifacts, such as halos, glare or shadowing, a major complaint of post capsulotomy patients.

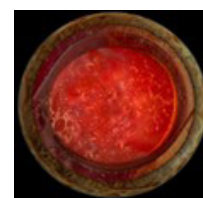

Actual view of the capsular opening using a Brasse WF Capsulotomy lens

For more information on this product, contact our customer service department.
Be sure to ask about our 30 day no risk trial!